

Operating conditions

- Do not operate the instrument in a moving vehicle
- Operating temperature 15–35°C for the instrument, 18–25°C for the test reagents
- Instrument storage at 2–35°C
- Reagent storage at 2–25°C

Contact and orders

Aidian Netherlands
info@aidian.nl / +31 317 351838
Orders: cs@aidian.nl
quikread.com

Aidian Belgium
info@aidian.be / +31 317 351838
Orders: cs@aidian.be
quikread.com

Product	Cat. no.
QuikRead go Instrument	OD-133893
QuikRead go easy CRP, including sample collectors, 50 tests	OD-153287
QuikRead go wrCRP+Hb, including capillaries and plungers, 50 tests	OD-146255
QuikRead go Bag, Torba TRM-42	OD-153902
Battery for QuikRead go Instrument	OD-136205

AIDIAN

aidian.nl / aidian.be

9194-02EN for Benelux, 08/2023. Not intended for the US market.

When minutes matter.

QuikRead go[®] is your lab on the move



Fast CRP test supports right treatment decision on the go!

QuikRead go[®] system is suitable for near-patient testing, for example in mobile units, home care units, pop-up testing sites, and ambulatory use.

Accurate and rapid CRP result aids evaluation of infection severity and supports fast decision making whether a patient should be sent to a hospital.



QuikRead go[®] easy CRP, used with the unique one-step sample collector, provides quantitative CRP result with minimal hands-on steps.



QuikRead go[®] wrCRP+Hb provides added value with combined quantitative wide measuring range CRP and hemoglobin results from a single sample.




Fast

- Quantitative CRP/CRP+Hb results in 2 minutes
- Small sample volume (10 µl) enables fast and easy testing
- Short start-up time of the instrument



Reliable

- Accurate results at POC – comparable to laboratory results¹
- Over 60,000 QuikRead go instruments in use worldwide

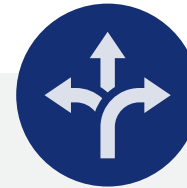

Flexible

- Sample is stable in the test buffer for 2 hours, either carry the instrument or the sample!
- The ready-to-use reagents can be stored at 2–25°C
- Suitable for near-patient testing in different use environments



Portable

- Long-lasting battery available
- Lightweight and portable instrument – carrier bag available
- Various connectivity solutions

Reference:

1. Instructions for use (QuikRead go[®] easy CRP / QuikRead go[®] wrCRP+Hb)